



BIOLOGICAL EVALUATION TEST REPORT

For

Tattoo Cartridge Needles / PMU&SMP Cartridge Needles

**Model: XXXXRL, XXXXRLT, XXXXRLMT, XXXXRLB, XXXXHRL, XXXXRS,
XXXXRS-1, XXXXRS-2, XXXXM1C, XXXXM1C-1, XXXXM1C-2, XXXXM1,
XXXXM1-1, XXXXM1-2, XXXXF, XXXXSF. (X=0,1,2.....9)**

Brand Name: EZ / EZ TATTOO / INKin / POPU

Report No.: ENC2308281GZ48E1

Date of Issue: Aug. 28, 2023


Prepared For

**Guixi City Jinlong Tattoo Equipment Manufactory
Economic Development Zone, Guixi City, Jiangxi Province , China**

Prepared By

**East Notice Certification Service Co., Ltd.
1/F, Haohui Commercial Building, Zhuji Street, Dongpu Town, Tianhe District,
Guangzhou City, China**

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

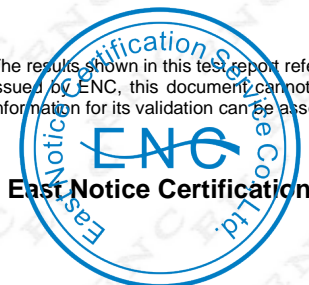

1/F, Haohui Commercial Building, Zhuji Street, Dongpu Town, Tianhe District, Guangzhou City
E-mail: enc@enc-lab.com Http:// www.enc-lab.com

VERIFICATION OF CONFORMITY

Product Designation:	Tattoo Cartridge Needles / PMU&SMP Cartridge Needles
Model:	XXXXRL, XXXXRLT, XXXXRLMT, XXXXRLB, XXXXHRL, XXXXRS, XXXXRS-1, XXXXRS-2, XXXXM1C, XXXXM1C-1, XXXXM1C-2, XXXXM1, XXXXM1-1, XXXXM1-2, XXXXF, XXXXSF. (X=0,1,2.....9)
Model Difference:	The series models have same materials as XXXXRL, except for the different appearance.
Brand Name:	EZ / EZ TATTOO / INKin / POPU
Applicant:	Guixi City Jinlong Tattoo Equipment Manufactory Economic Development Zone, Guixi City, Jiangxi Province , China
Manufacturer:	Guixi City Jinlong Tattoo Equipment Manufactory Economic Development Zone, Guixi City, Jiangxi Province , China
Type of Test:	DELAYED-TYPE HYPERSENSIVITY TEST
Technical Standards:	EN ISO 10993-1:2020, EN ISO 10993-10:2021, EN ISO 10993-12:2021
File Number:	ENC2308281GZ48E1
Date of test:	Jul. 27, 2023 – Aug. 28, 2023
Deviation:	None
Condition of Test Sample:	Normal

The above equipment was tested by East Notice Certification Service Co., Ltd. For non sterile medical device Class I according Annex VII of the council directive 93/42/EEC and the Technical Standards mentioned above. Should any objections to the test reports occurred, should submit it to the Company within ten days since the issuing of the report, Fail to accept.


The test results of this report relate only to the tested Sample identified in this report.

Checked By Yemig
Yemig Aug. 28, 2023

Authorized By Ray Zhou
Ray Zhou Aug. 28, 2023




The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

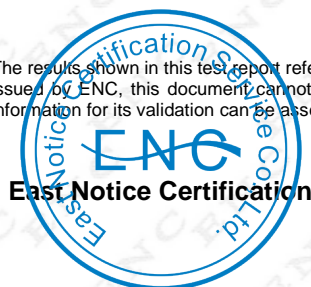

INDEX

SUMMARY	3
INTRODUCTION	4
BIBLIOGRAPHY	4
FILING	4
PROCEDURES	5
TEST SUBSTANCE DESCRIPTION	5
ANALYSED SAMPLE	5
DELAYED-TYPE HYPERSENSIVITY TEST	6
EXPERIMENTAL PROCEDURE	6
INTERPRETATION OF RESULTS	10
RESULTS	10
CONCLUSIONS	10
APPENDICES	11

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

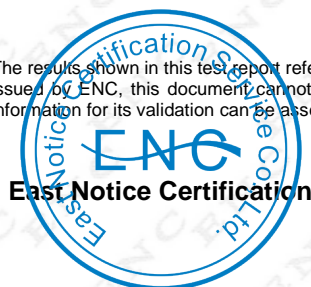

SUMMARY

A toxicological study was performed to evaluate the biocompatibility of the test substance Tattoo Cartridge Needles / PMU&SMP Cartridge Needles, at this purpose the following test was carried out:

- test standards: EN ISO 10993-1:2020, EN ISO 10993-10:2021, EN ISO 10993-12:2021

The analytical test was accomplished on the materials of Needle and Pc Plastic parts which constitute the device and are in contact with the human skin:

Two eluates of the test substance were prepared both in vegetable oil and in physiological solution in order to perform the delayed-type hypersensitivity test.

In static condition the eluates of the test material were performed by immersing the test material in both physiological solution and vegetable oil in order to reach weight/volume ratios and surface/volume ratios of:

- 0.2 g /ml for the PC Plastic parts
- 0.2 g/ml for the Needle

The test sample was then incubated for 72 hours at $37^{\circ}\text{C} \pm 1^{\circ}\text{C}$, after this period, has been done a pool of eluates.

For each elution 15 guinea pigs were used, 10 treated with the eluate of the test substance and 5 using as control treated with extraction liquid only.

The skin sensitisation test had 2 phases, induction phase and challenge phase.

During the induction phase the group of 10 treated guinea pigs were treated with 3 double intradermal injections as follows:

- 1st Freud Complete Adjuvant in distilled water (ratio 1:1)
- 2nd Elution of the test substance
- 3rd Elution of the test substance and FCA (ratio 1:1)

The control animals received the same pairs of injections, but in the 2nd injection only extraction liquid was administered (physiological solution or vegetable oil).

In the third injection, extraction liquid + FCA (ratio 1:1) was used.

After 3 days from the beginning of treatment on the all animals, a topical application, with slight massage, of 0.5 ml of Sodium Lauril Sulfatum 10% was performed.

After 4 days from the intradermal injections, the test substance was applied (at a dose of 0.5 ml/animal). The application lasted 48 hours.

The same treatment was used on control guinea pigs using only extraction liquid.

After 7 days from the beginning of treatment the challenge phase was performed by applying 0.5 ml of the eluate on the left side and 0.5 ml of the eluates on the right side. The

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.



bandage was left on for 24 hours. 24, 48 and 72 hours after the beginning of this phase, the tested and the control animals were observed.

Neither edema nor erythema were observed in the animals treated with the test substance eluate prepared in physiological solution. and in vegetable oil.

No abnormalities were observed in the animals used as control.

On the basis of the results obtained for all tested components, interpreted according to EN ISO 10993-10:2021 the test substance Tattoo Cartridge Needles / PMU&SMP Cartridge Needles must be considered **NOT SENSITIZING**.

INTRODUCTION

This study has been carried out on behalf of Needle and PC Plastic parts on the product.

The study was performed at the Test Facility ENC-lab of East Notice Certification Service Co., Ltd.

TEST	START	END	RESEARCHER
delayed-type hypersensitivity test	Jul. 27, 2023	Aug. 28, 2023	Sam Liu

BIBLIOGRAPHY

1. EN ISO 10993-10:2021

Biological evaluation of medical devices

Part 10: Tests for irritation and delayed-type hypersensitivity


FILING

The study program, all raw data and a copy of the final report are filed in the archives of ENC-lab. for ten years after the issuing of the final report.

No retained sample will be kept.

At the end of the conservation period, the Sponsor may request an extension of the conservation of all or part of the substances for a further period, or their restitution. A suitable agreement shall be drafted in this case.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

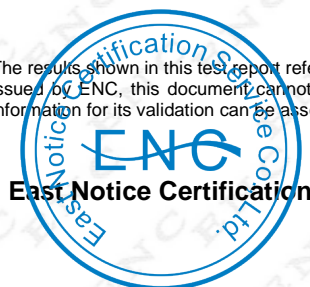

PROCEDURES

All procedures used during this study are recorded in the ENC-lab Procedures Manual.

TEST SUBSTANCE DESCRIPTION

The test substance is a device consisting of Needle and PC Plastic parts intended to human use in contact with the skin.


Name: Tattoo Cartridge Needles / PMU&SMP Cartridge Needles

ANALYSED SAMPLE

The analysed sample, representative of the test substance, is identified by the following numbers:

Name:	Tattoo Cartridge Needles / PMU&SMP Cartridge Needles
Test Model:	XXXXRL
Date of production:	Jun. 22, 2023
Shelf life:	2 years
Receiving date:	Jul. 27, 2023
Registration number:	N/A

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

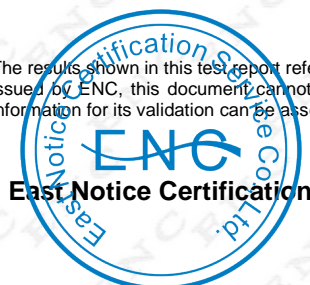

DELAYED-TYPE HYPERSENSIVITY TEST

EXPERIMENTAL PROCEDURE

1. TEST METHOD

1.1 Characterisation

Species: Albino guinea pigs
Strain: Hartley
N.: 30
Weight: 300 - 400 g at the arrival at the Centre
Sex: female
Supplier: Guangzhou Juyuan Breeding Farm

1.2 Caging

The animals were caged, in groups of ten, in transparent polycarbonate cages (dimensions: 590x385x200h mm).

The housing room was lighted with fluorescent lamps 12 hours for day.

Room temperature and humidity were regulated by a conditioning plant and were monitored daily.

Recordings of the housing conditions are being retained in ENC-lab files.

1.3 Cleaning and disinfection

The cages and the housing room were cleaned before the animals were accommodated, then disinfected periodically.

1.4 Feeding

Animals have been fed with standard pellet complete diet supplied by the authorized breeder Susan.

1.5 Watering

Filtered tap water from local network was supplied ad libitum from an automatic watering system.


1.6 Quarantine

Before allocation to the study, the animals were kept in quarantine for one week.

During this period they were observed daily.

At the end of the quarantine period the animals were carefully examined in order to evaluate their suitability for the study.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

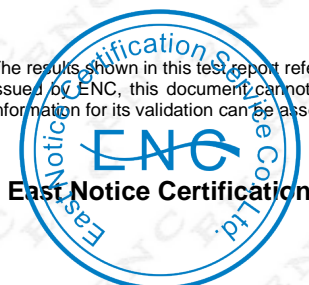

2. TEST SAMPLE PREPARATION

The eluates of the test substance were prepared in static conditions by immersing:

- PC Plastic parts (22g) into 110 ml of both eluants in order to reach a weight/volume ratio of 0.2 g/ml.
- Needle (60g) into 300 ml of both eluants in order to reach a weight/volume ratio of 0.2 g/ml.

The test samples were then incubated for 72 hours at 37⁰C ±1⁰C, after this period, two pools of eluates have been done.

3. EXPERIMENTAL DESIGN

Experimental design consisted of two groups (treated) of 10 animals treated with extract in vegetable oil and in physiological solution of test substance (group 1-2) and two group (control) of 5 control animals treated with only vegetable oil or physiological solution (group 3-4). The animals were divided in groups as follows:

GROUP	INDUCTION		CHALLENGE TOPIC APPLICATION
	Intradermal injection	Topic application	
1	1. Extract in physiological solution 2. Extract in physiological solution + FCA 3. FCA	Extract in physiological solution	Right side: Extract in physiological solution
			Left side: Physiological solution
2	1.Extract in freshly refined vegetable oil 2.Extract in freshly refined vegetable oil+ FCA 3.FCA	Extract in vegetable oil	Right side: Extract in vegetable oil
			Left side: Vegetable oil
3	1. Physiological solution 2. Physiological solution + FCA 3. FCA	Physiological solution	Right side: Extract in physiological solution
			Left side: Physiological solution
4	1. Freshly refined vegetable oil 2. Freshly refined vegetable oil + FCA 3. FCA	Vegetable oil	Right side: Extract in vegetable oil
			Left side: Vegetable oil

The animals allocated to the study were selected randomly from those suitable, available at that time.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.



4. TREATMENT

4.1 Skin preparation

24 hours before testing, fur was removed by shaving a 50 cm² wide area on the back of the animals.

4.2 Administration

The test consisted of an induction phase and a challenge phase.

Induction phase

Day 0 - treated group

Three pairs of 0,1 ml intradermal injections were made in the intrascapular region of each animal, on each side of the midline, according to the following scheme:

- 1) FCA in distilled water (ratio 1:1)
- 2) Elution of test substance
- 3) Elution of test substance + FCA (ratio 1:1)

Day 0 - control group

Three pairs of 0,1 ml intradermal injections were made in the intrascapular region of each animal, on each side of the midline. The content was:

- 1) FCA in distilled water
- 2) Extraction liquid
- 3) Extraction liquid + FCA (ratio 1:1)

Day 3 - treated group and control group

After 6 days the beginning of treatment on the all animals a topical application, with slight massage of 0,5 ml of Sodium Lauril Solfatum 10%, was made.

Day 4 - treated group

Seven days after the intradermal injections had been made, 0,5 ml of the elution of the test substance were applied to each animal and held in place with an occlusive patch. The application was made on area caudally to the area of injection.

The dressing was left in place for 48 hours.

Day 4 - control group

The same treatment was performed on the control group, using vegetable oil and physiological solution instead of the test substance.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.



Challenge

Day 7 - treated and control groups

An occlusive patch with 0,5 ml of eluate of the test substance was applied to the right side and to the solvent to the left side.

The dressing was left in place for 24 hours.


OBSERVATIONS

On the 3rd day (24 hours after removal the patch), and the 4th day (48 hours after removal the patch) and the 5th day (72 hours after removal the patch) of tests all the animals treated and controlled were evaluated for a skin reaction.

The intensity of erythema and/or edema were evaluated according to the following scale:

Reaction	Grade
<u>Erythema</u>	
No erythema	0
Slight erythema	1
Well defined erythema	2
Moderate erythema	3
Severe erythema to slight eschar formation	4
<u>Edema</u>	
No edema	0
Slight edema	1
Well defined edema	2
Moderate edema	3
Severe edema	4

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

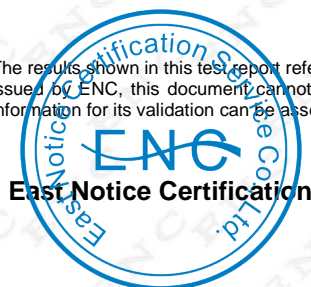

INTERPRETATION OF RESULTS

Both frequency and intensity of response were evaluated.

In case of positive reaction only in treated animals, the frequency of sensitization was considered, without taking into account the intensity of the response.

RESULTS

Neither edema nor erythema were observed in the animals treated with the test substance eluates prepared in physiological solution and in vegetable oil.

No abnormalities were observed in the animals used as control.

% sensitising guinea pigs treated with extract in physiological solution: 0%


% sensitised guinea pigs treated with extract in vegetable oil: 0%

The data concerning every single animal are reported in appendices 1 and 2.

CONCLUSIONS

On the basis of the results obtained for all tested components, interpreted according to EN ISO 10993-10:2021 the test substance Tattoo Cartridge Needles / PMU&SMP Cartridge Needles must be considered **NOT SENSITIZING**.

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

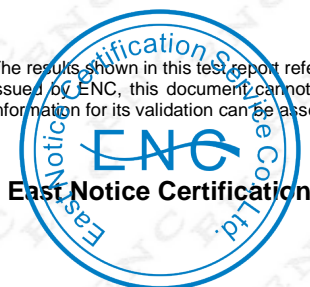

APPENDICES

Appendix n.1: Skin reactions in treated animals with eluate

ANIMAL N.	TIME AFTER CHALLENGE APPLICATION					
	physiological solution					
	48 hours		72 hours		96 hours	
	Erythema	Edema	Erythema	Edema	Erythema	Edema
1	0	0	0	0	0	0
2	0	0	0	0	0	0
3	0	0	0	0	0	0
4	0	0	0	0	0	0
5	0	0	0	0	0	0
6	0	0	0	0	0	0
7	0	0	0	0	0	0
8	0	0	0	0	0	0
9	0	0	0	0	0	0
10	0	0	0	0	0	0

ANIMAL N.	TIME AFTER CHALLENGE APPLICATION					
	vegetable oil					
	48 hours		72 hours		96 hours	
	Erythema	Edema	Erythema	Edema	Erythema	Edema
1	0	0	0	0	0	0
2	0	0	0	0	0	0
3	0	0	0	0	0	0
4	0	0	0	0	0	0
5	0	0	0	0	0	0
6	0	0	0	0	0	0
7	0	0	0	0	0	0
8	0	0	0	0	0	0
9	0	0	0	0	0	0
10	0	0	0	0	0	0


Erythema

0= No erythema

Edema

0= No edema

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

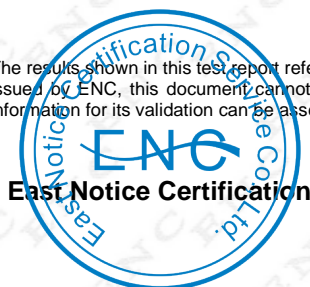

Appendix n.2: Skin reactions in control animals treated with vegetable oil and physiological solution

ANIMAL N.	TIME AFTER CHALLENGE APPLICATION					
	physiological solution					
	48 hours		72 hours		96 hours	
	Erythema	Edema	Erythema	Edema	Erythema	Edema
1	0	0	0	0	0	0
2	0	0	0	0	0	0
3	0	0	0	0	0	0
4	0	0	0	0	0	0
5	0	0	0	0	0	0

ANIMAL N.	TIME AFTER CHALLENGE APPLICATION					
	vegetable oil					
	48 hours		72 hours		96 hours	
	Erythema	Edema	Erythema	Edema	Erythema	Edema
1	0	0	0	0	0	0
2	0	0	0	0	0	0
3	0	0	0	0	0	0
4	0	0	0	0	0	0
5	0	0	0	0	0	0


Erythema

0= No erythema

Edema

0= No edema

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

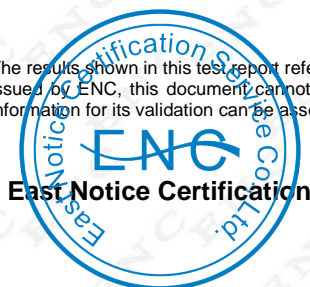

Appendix n.3: Photographs of sample



---- END OF REPORT ----

The results shown in this test report refer only to the sample(s) tested unless otherwise stated and the sample(s) are retained for 30 days only. The document is issued by ENC, this document cannot be reproduced except in full with our prior written permission. The document is available on request and the brief information for its validation can be assessable and confirmed at <http://www.enc-lab.com>.

